

DAC[®]

DEFENSIVE
ANTIBACTERIAL
COATING

ADLER[®]
ORTHO



contents

1. Implant risk
2. Bacteria mechanism & biofilm stages
3. Defensive Antibacterial Coating - DAC®
4. DAC® difference
5. Relative risks
6. DAC® features and benefits
7. Prophylactic infection prevention
8. Mixing DAC®
9. Any implant DAC®
10. DAC® and antibiotics
11. DAC® evidence and references

implant risk

Biofilms can form on implants in the first few post operative hours and can determine the ultimate outcome of the procedure.

Introducing foreign materials, like orthopaedic implants into the body carries the risk of infection due to possible biofilm formation. Biofilms are a particular concern because they resist antibiotics and host defense systems.

Despite advancements in orthopaedic devices and surgery, managing infection remains challenging. Infections invariably lead to prolonged treatment, unstable fixation, early revision, implant removal, functional loss, and sometimes amputation. These consequences strain hospitals financially and greatly impact patients' well-being.

Although global infection rates are low in the primary population, peri-prosthetic joint infection (PJI) and fracture-related infection (FRI) can occur more frequently in patients with increased risk factors. ^{1, 2}

While the factors contributing to implant infections may vary, focusing on prevention of biofilm formation can mitigate the environmental impact, resource burdens (such as bed space) and overall direct end costs. ³

The impact on patients, their carers and families as well as the clinical team should not be underestimated.

bacteria mechanism

There are more bacteria cells than human cells in our bodies ⁴. Bacteria typically exist freely in a planktonic state within the body. However, when they encounter an implant surface, they can attach and form communities. Within hours, these bacteria enhance their survival by encasing themselves in exopolysaccharides, forming a biofilm on the **entire implant surface, including at the bone-implant interface.**


The formation of this biofilm can negatively impact the fixation and overall stability of the implant, potentially leading to the need for revision surgery. Once established, biofilms are notoriously resistant to removal, withstanding shear forces and antibiotic treatments, which makes eradication challenging.

race for the surface

The concept of “race for the surface,” highlights the competition between microbial adhesion and tissue integration on biomaterial implants. Typically, the crucial phase occurs within the first 72 hours. If tissue cells win, the implant is covered by tissue and remains resistant to bacterial colonization. If bacteria win, a biofilm forms, hampering tissue cell function and increasing infection risk.

biofilm stages

The biofilm formation process is integral to the five stages of infection: ⁵

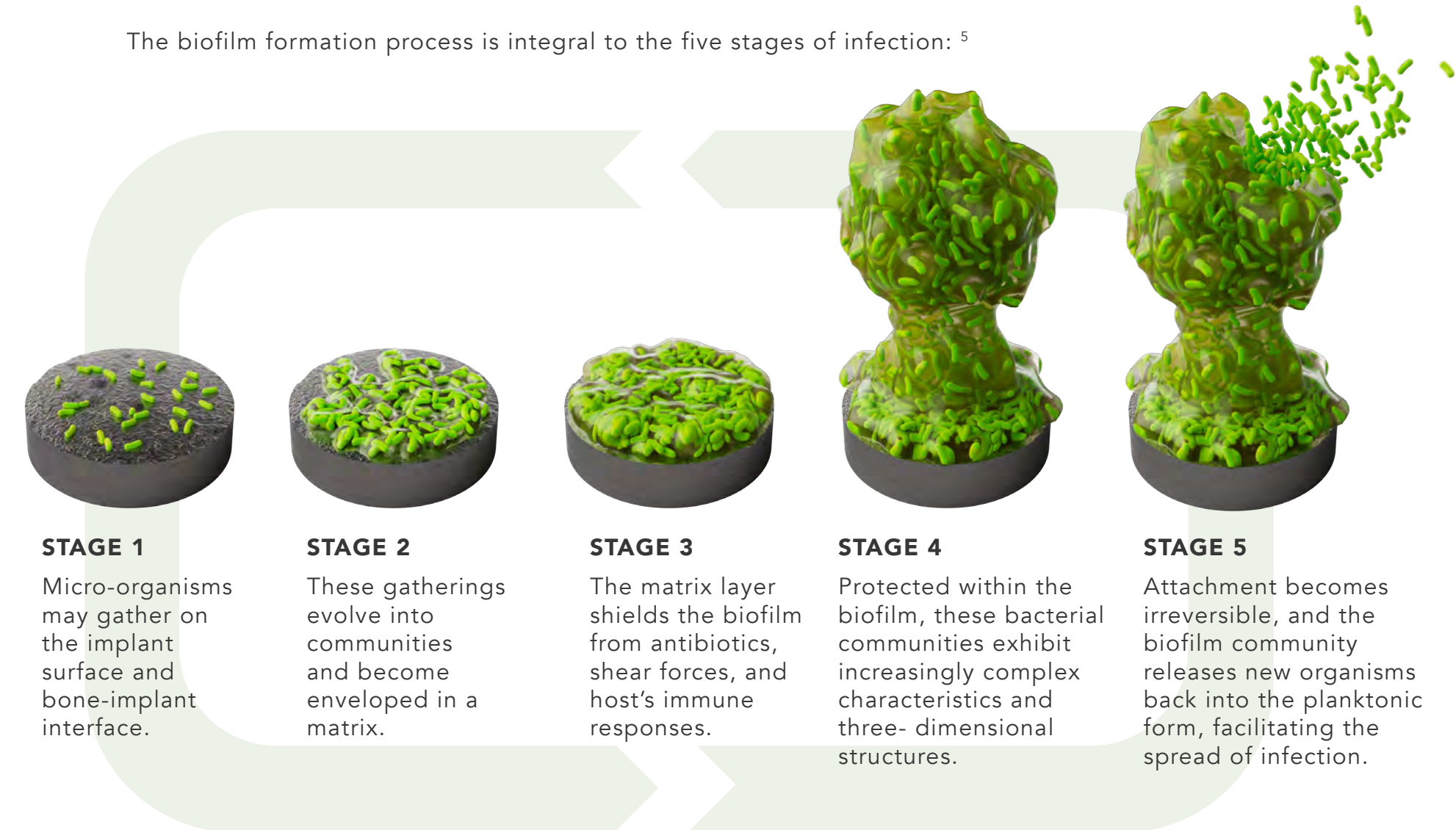

DAC[®] defensive antibacterial coating

Defensive Antibacterial Coating (DAC[®]) is a smart hydrogel with hydrophilic properties, which, when spread onto the implant surface, acts as a physical barrier to prevent biofilm formation, defending against potential post-surgical infections, and optimising patient outcomes.

DAC[®] is a powdered medical device comprised of two bioresorbable polymers, designed for use in orthopaedic surgeries. DAC[®] forms a temporary hydrogel that adheres seamlessly to implant surfaces. This innovative hydrogel serves as a prophylactic coating, which can be easily applied both before and after the implantation of orthopaedic devices.



HYALURONIC ACID (HA)

Hyaluronic acid (HA) is a naturally occurring substance found in the human body. Highly biocompatible, surfaces coated with HA show minimal bacterial biofilm growth, making it renowned for its exceptional moisture retention capabilities.

POLY-LACTIC ACID (PLA)

Poly-lactic acid (PLA) is a safe, biodegradable and bio-absorbable synthetic polymer obtained from renewable sources (corn or other cereals). It is commonly used in various applications, including medical devices. PLA decomposes into harmless byproducts in the environment, contributing to its eco-friendliness.

Surgeons may choose to load antibiotics to DAC[®].

DAC[®] difference


DAC[®] effectively prevents biofilm formation on the implant surface, including the bone-implant interface, during the critical first 72 hours, reducing infection risks in orthopaedic procedures. It allows surgeons to coat any implant with a hydrogel, preventing bacterial adhesion and enabling the direct application of antibiotics to the site.

DAC[®] consistently demonstrates effectiveness across various challenging surgical procedures, facilitating localized antibiotic release to target infection-prone areas.

Clinical assessments, such as those for periprosthetic joint infections in total hip revision arthroplasty ¹³, support its efficacy. Using antibiotic-loaded DAC[®] with cementless prostheses, alongside systemic prophylaxis, reduces early postoperative infections. Furthermore, a comprehensive review of clinical data ⁶ highlights DAC[®]'s efficacy, with notable success in high-risk vertebral surgeries.


These findings underscore DAC[®]'s potential in enhancing infection control and improving surgical outcomes.

without DAC[®] coating applied

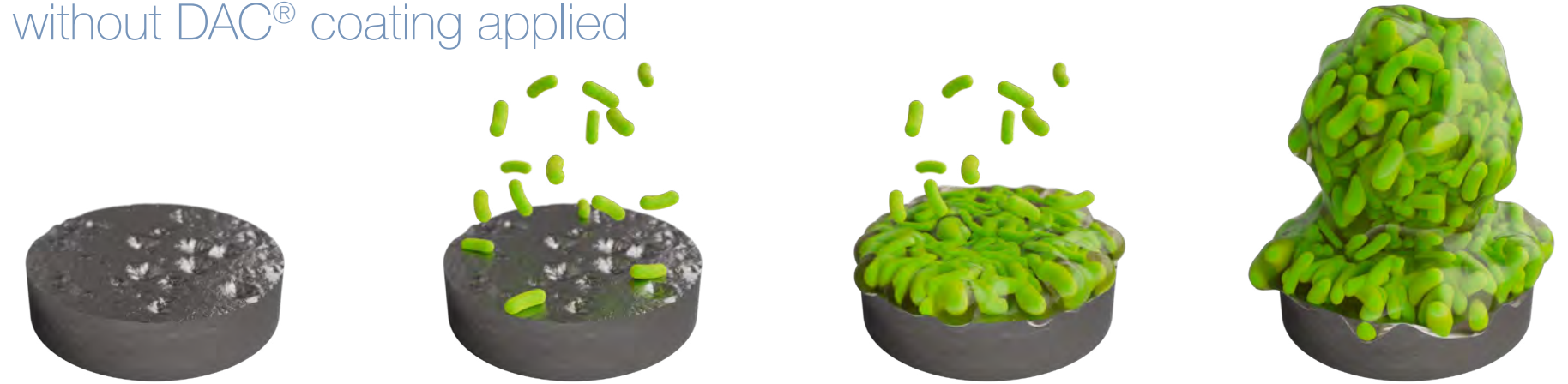

Implants are susceptible to pathogens adhering to their surface. In the absence of a response to the organism, bacterial colonies initiate biofilm production.

DAC[®] coating applied

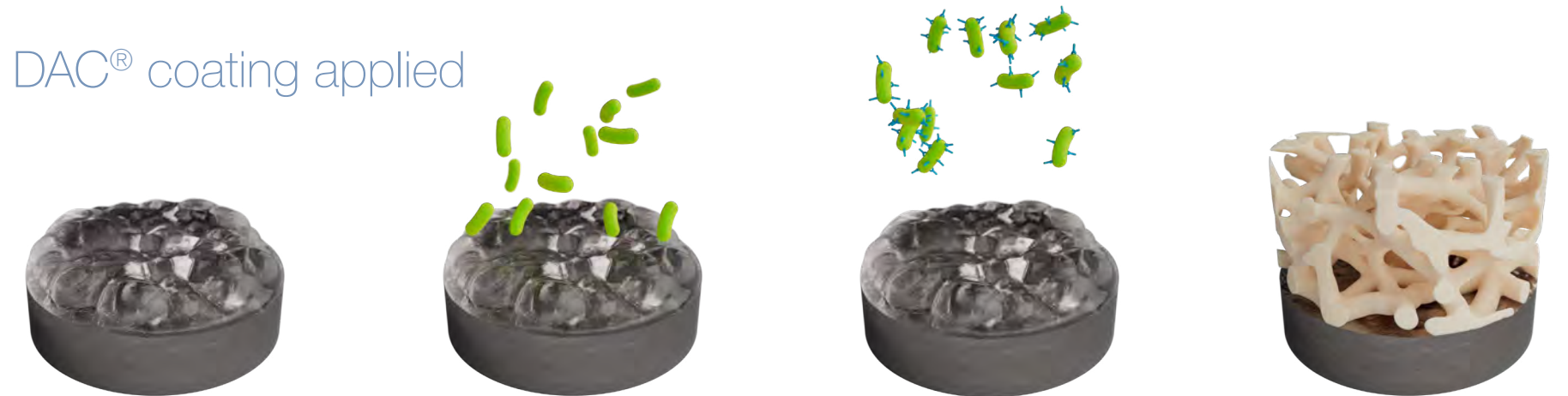

The DAC[®] hydrogel layer interrupts bacterial adhesion to the implant surface, due to its hydrophilic properties. The immune system and administered systemic antibiotics can effectively target planktonic micro-organisms.

relative risks

In any surgical setting, infection prevention is paramount, even among healthy individuals. While the relative risks in primary settings are generally low, they can quickly escalate with the presence of one or more risk factors. These factors, whether in primary, revision joint replacement, trauma, or oncology cases, increase the patient's susceptibility to infection.

The true value of DAC® becomes particularly evident in cases involving higher-risk patients, complex surgeries, or when additional factors need consideration. In these scenarios, the benefits of DAC® can significantly outweigh its costs, especially given the increasing prevalence of such cases.¹⁴




other factors:

An observational cohort study in the UK analysed 623,253 primary hip procedures, focusing on revisions due to periprosthetic joint infections (PJI). Data from the National Joint Registry, linked to Hospital Episode Statistics, identified several risk factors for PJI after primary hip replacement:

- Elevated BMI
- Comorbidities such as diabetes
- Previous septic arthritis or infection
- Previous femoral neck fracture

These factors are associated with an increased risk of PJI, highlighting the need for targeted infection control measures in these patient groups.¹⁵

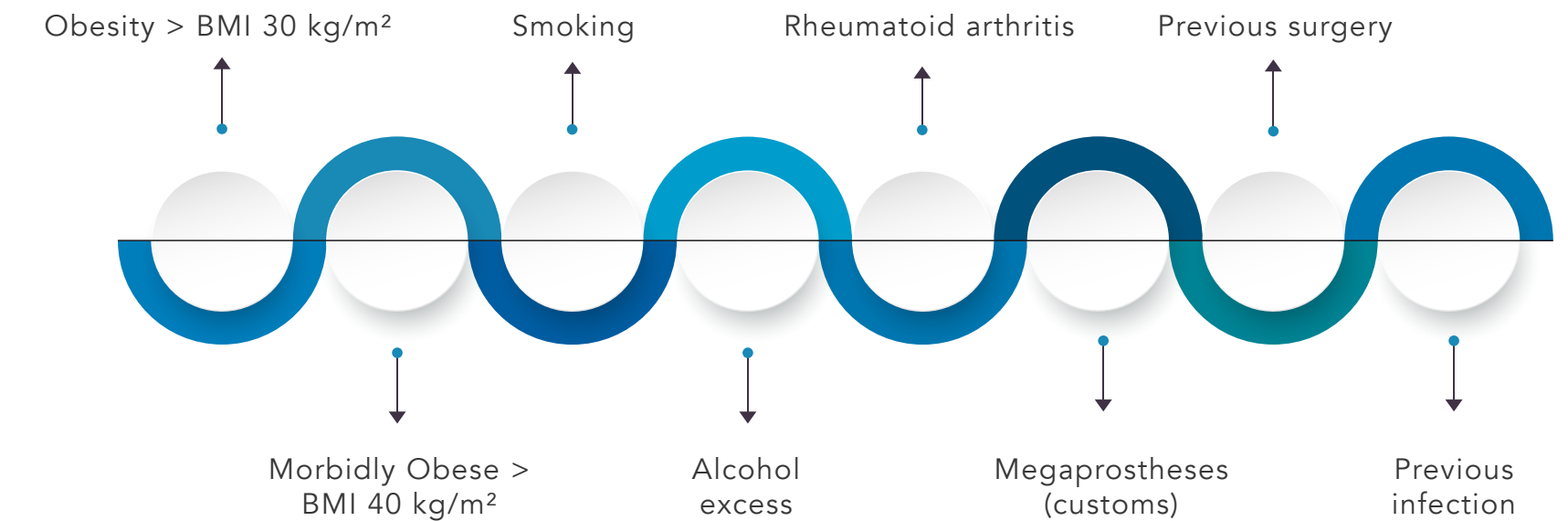

Furthermore, a large-scale analysis¹⁶ involving over 2 million patients identified higher BMI, particularly obesity (BMI above 30 kg/m²), as a significant risk factor for both deep and superficial infections. Morbidly obese individuals (BMI over 40 kg/m²) exhibited an even higher susceptibility. Additional risk factors included:

- High ASA grade
- Trauma site surgery
- Certain ethnic groups


- Comorbidities
- Previous infections
- Multiple site operations
- Use of mega-protheses
- Smoking

These factors collectively increase the risk of infections, emphasizing the importance of tailored infection control strategies in orthopaedic surgery.

DAC[®] features and benefits

Numerous investigations, surgeon congresses, and initiatives have sought to understand and alleviate the pressures and financial challenges of orthopaedics infections, with varying degrees of success.

However, some efforts have resulted in direct and indirect issues such as wound leakage, third body wear, mechanical limitations, off-label product use, thermal damage to antibiotics, and migration from the intended site, leading to complications during re-revision.

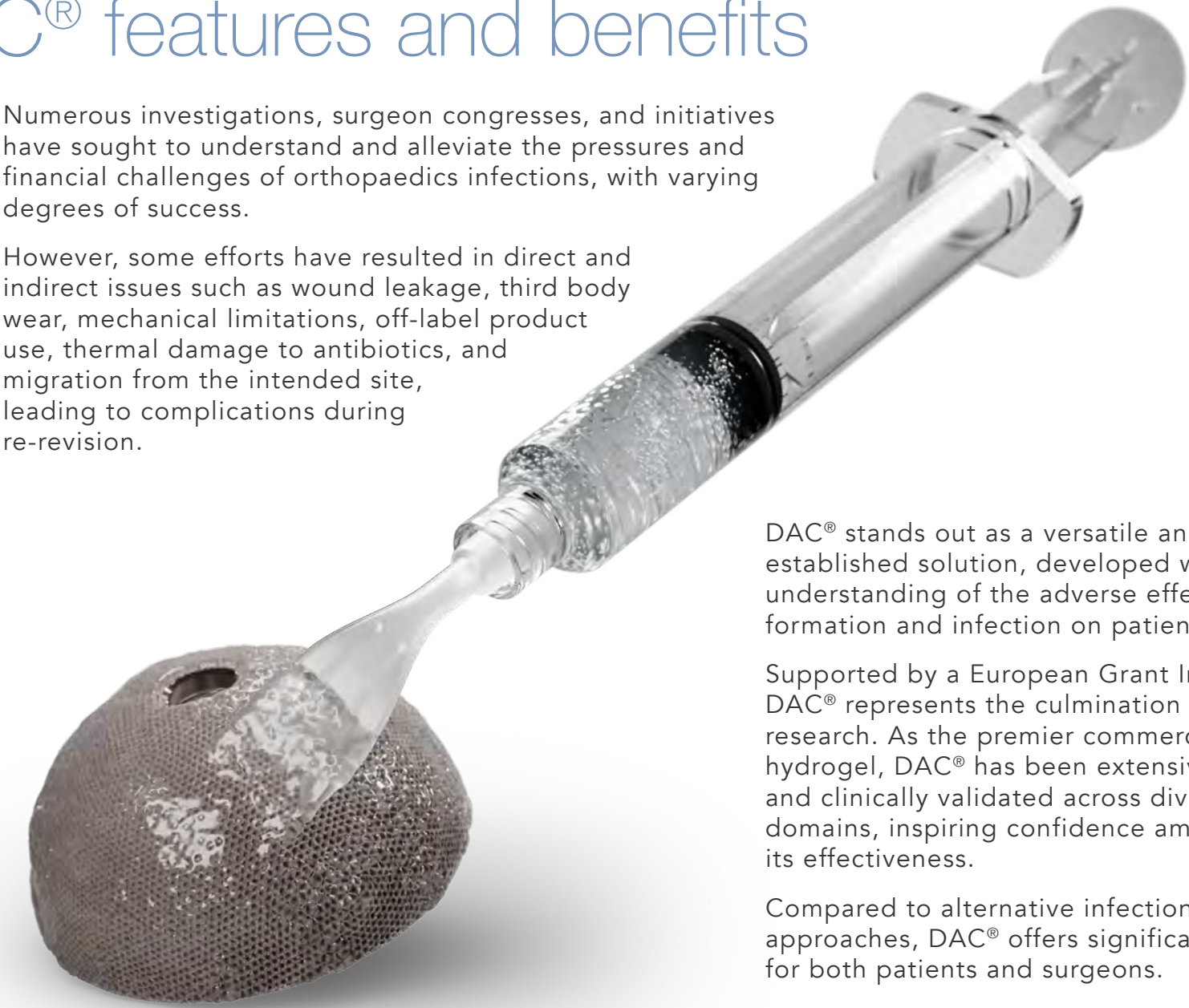

DAC[®] stands out as a versatile and well-established solution, developed with a deep understanding of the adverse effects of biofilm formation and infection on patient outcomes.

Supported by a European Grant Initiative, DAC[®] represents the culmination of extensive research. As the premier commercially available hydrogel, DAC[®] has been extensively published and clinically validated across diverse anatomical domains, inspiring confidence among clinicians in its effectiveness.

Compared to alternative infection prevention approaches, DAC[®] offers significant advantages for both patients and surgeons.

-  Indicated for prevention of infection ✓
-  Clinically proven ✓
-  Articulating joints ✓
-  Can be used with trauma plates and nailing ✓
-  Cemented and cementless devices ✓
-  Easy to mix ✓
-  Simple to apply ✓
-  Can be used with antibiotics ✓
-  Active for 72 hours ✓
-  Suitable for trauma, revision, oncology implants ✓

DAC[®] prophylactic infection prevention

DAC[®] has proven its capacity to improve overall outcomes and healthcare economics with substantial clinical evidence spanning the past decade. Its impact is particularly notable in reducing infection-related costs, which can soar for revision cases. This cost reduction can be achieved through DAC[®]'s effectiveness, demonstrating efficacy even at critical interfaces like the bone implant, in articulating joint replacements and with cementless implants in long bones.


DAC[®] presents a reliable infection prevention strategy across a spectrum of surgical scenarios, including trauma, revision, elective, and oncology procedures. Surgeons rely on DAC[®] with confidence, as its safety and efficacy have been clearly demonstrated through post market surveillance.

Beyond infection prevention, DAC[®] completely disaggregates and does not hinder implant osseointegration, and ensures implant adherence, rendering it suitable for comprehensive coating of entire devices.

Notably, DAC[®] is specifically designed for infection prevention, distinguishing it from other strategies. This specialisation further enhances its effectiveness in safeguarding patient health and minimising the burden of post-operative infections.

DAC[®] is effective in various surgical scenarios, including:

- Trauma
- Elective procedures
- Oncology surgeries
- Cemented or cementless procedures
- Megaprotheses
- DAIR cases
- Nailing & Plates
- Articulating joint replacements
- Prevention of fracture-related infections in both primary and revision surgeries

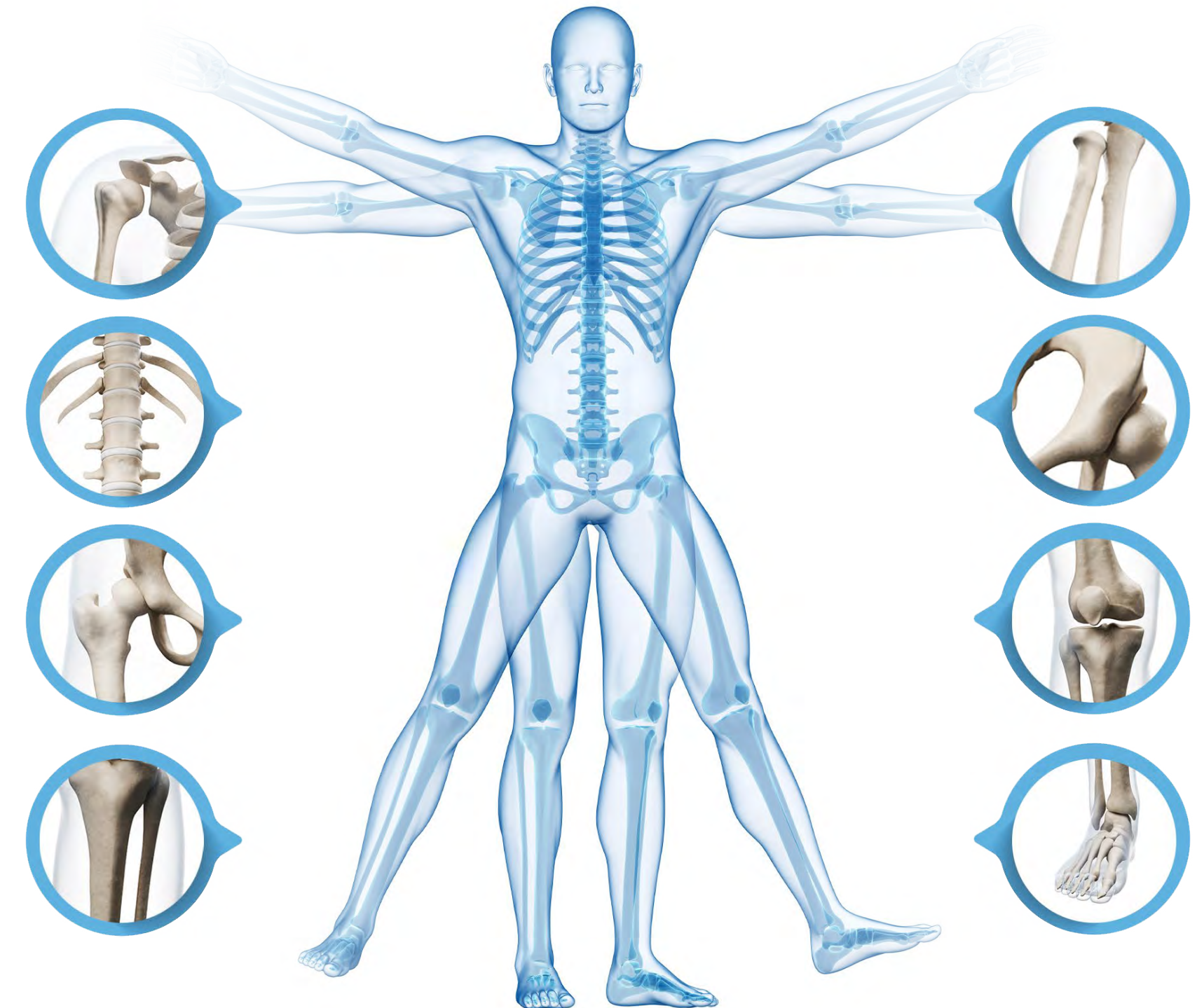

mixing DAC®

1



Before use, de-compact the powder as necessary by carefully and slightly retracting the syringe plunger position and digitally flicking or tapping the syringe until the powder is loosened and moving freely.

2



Connect the backstop to the DAC® powder syringe. Ensure the BD is pointing up towards the plunger.

3



Draw up 5ml of sterile water (or 5ml of antibiotic solution).

4




With the syringe in an upright vertical position and plunger towards table, unscrew the cap from the DAC® syringe. Screw on the luer-lock connector.

5




Connect the solution syringe and place the fluid in the top and powder at the bottom.

6

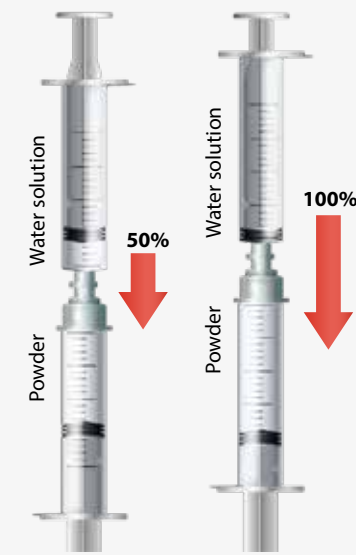

Press 50% of the solution firmly and quickly into the DAC® powder. Press the remaining 50% of fluid slowly into the powder maintaining a vertical position of the syringe throughout this motion.

7

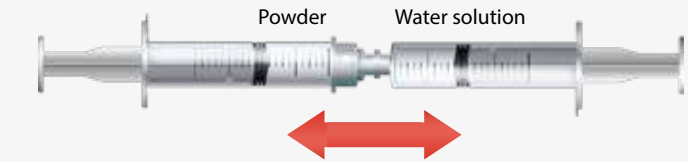

Place the syringe into a horizontal position and move the solution part way between the two for the first few passes. Once the mix starts to take place then fully depress each syringe until a homogenous mix is achieved.

8



Mix is complete once all powder is dissolved. Ensure the hydrogel is transferred completely into one syringe. Disconnect syringes. Screw on the applicator and leave for 5 to 10 minutes minimum before application. Once DAC® is prepared it is safe to rest for up to 6 hours before use.

9



The implant surface is coated with DAC® preventing bacterial adhesion and colonisation of the implant.

any implant DAC[®]

DAC[®] should be kept cool and stored in a fridge between 2°C and 8°C to ensure it is easier to mix and the HA is stable.

Unmixed DAC[®] should not be used if the ambient temperature has reached 25°C since the HA can become unstable.

Surgical Procedure:

Mix DAC[®] at the very start of any procedure before the skin incision. This ensures the mix is ready when required and once mixed into a stable compound the temperature is no longer of concern.

In any event DAC[®] should be mixed 15 minutes before it is to be applied to an implant.

Once mixed and hydrated DAC[®] can safely be used for up to 6 hours.

To prepare the DAC[®] :

DAC[®] Kit
Sterile water
Small sterile pot

Antibiotic if required
Syringe needle to draw up antibiotic



laboratory studies

Laboratory studies have demonstrated high rate of infection eradication following cementless one-stage revision hip arthroplasty using antibacterial hydrogel coating ¹⁷. The DAC[®] coating layer remains attached to the device in the case of a cementless stem ¹⁸ and DAC[®] does not impact the bony ingrowth of cementless implants. Moreover, the product can be applied between the implant bone interface safely.

DAC[®] demonstrates remarkable efficacy in numerous indications including cemented and cementless arthroplasty, including those implants which are already in situ during revision surgery. Do not use DAC[®] on an infected area not surgically cleaned.

Recommended techniques include:

For cementless devices: Coat all areas within bone or interfacing with bone. Complete any pulse lavage stages. Coat any remaining areas outside the bone. For these products the ingrowth has been shown to be unaffected. The hydrogel remains on the device during insertion due to its adhesive properties.

For cemented devices: Cement the devices into place, perform final pulse lavage, and coat all visible areas.

For products requiring assembly: Fully assemble any morse tapers before coating the device. Under no circumstances should the hydrogel be applied to the inside of a hip head.

Trauma plates: Coat the backside of the plate, affix to the bone, perform pulse lavage, and coat the remainder of the device.

The coating should be uniformly and consistently applied in a thin layer to ensure comprehensive implant protection.

Care should be taken in application of DAC[®] for operations involving the spine, which should only be applied to coat vertebral stabilisation devices, including screws, rods and cages.

(Please refer to DAC[®] IFU for further information.)

DAC® and antibiotics

DAC003000

Single Kit comprises:

- 1 sterile DAC® syringe containing 300 mg of dry powder
- 1 complete set of sterile DAC® components (connector, back-stop, spreader)
- 1 empty 10ml graduated syringe
Designed to prepare 5 ml of DAC® Hydrogel when reconstituted.

DAC003002

Double Kit comprises:

- 2 sterile DAC® syringes containing 300 mg of dry powder each
- 2 complete sets of sterile DAC® components (connector, back-stop, spreader)
- 2 empty 10 ml graduated syringes
Intended to prepare 10 ml of DAC® Hydrogel when reconstituted.


Storage:

DAC® products should be stored in a refrigerator at temperatures between 2°C and 8°C.

Do not freeze. Refer to the full Instructions for Use (IFU) for comprehensive instructions.

Clinical Usage:

DAC® has been utilised effectively in clinical settings both with and without antibiotics.

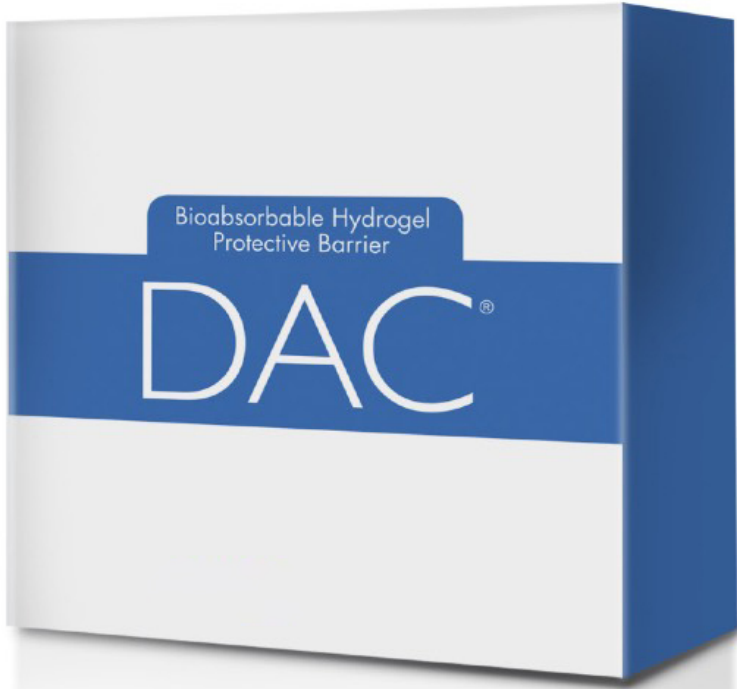

The following chart serves as a guide for surgeons who opt to use DAC® in combination with antibiotics, in regions where regulatory clearance permits.

Guide	DAC®	DAC® Powder	Vancomycin	Sterile water	Draw up Total	Resultant antibiotics per DAC®
Primary / trauma plate	DAC003000 - 5 ml	300 mg	1 g	20 ml	5 ml	250 mg / 5 ml DAC®
Primary / trauma plate	DAC003000 - 5 ml	300 mg	500 mg	10 ml	5 ml	250 mg / 5 ml DAC®
Revision / trauma nail	DAC03002 - 10 ml	600 mg	1 g	20 ml	10 ml	500 mg / 10 ml DAC®
Revision / trauma nail	DAC03002 - 10 ml	600 mg	500 mg	10 ml	10 ml	500 mg / 10 ml DAC®
For megaprosthes / larger implants additional DAC® may be required						

Antibiotic to be used	Antibiotic quantity	Ampoules (as applicable)	DAC® quantity	Volume of sterile water	Draw up volume	Total Antibiotics / volume
Vancomycin	1 g	-	300 mg - 1 syringe	20 ml	5 ml per syringe	250 mg / 5 ml DAC®
Vancomycin	500 mg	-	300 mg - 1 syringe	10 ml	5 ml per syringe	250 mg / 5 ml DAC®
Vancomycin	1 g	-	600 mg - 2 syringes	20 ml	5 ml per syringe = 10ml total	500 mg / 10 ml DAC®
Vancomycin	500 mg	-	600 mg - 2 syringes	10 ml	5 ml per syringe = 10ml total	500 mg / 10 ml DAC®
Gentamycin	80 mg/ 2 ml	2	300 mg - 1 syringe	1 ml	5 ml per syringe	160 mg/ 5 ml DAC®
Gentamycin	80 mg/ 2 ml	5	600 mg - 2 syringes	-	5 ml per syringe	200 mg / 5 ml DAC®
If the surgeon considers using other antibiotics if appropriate, then the final concentration should be not less than 2% and not more than 5%.						
Antibiotics are not supplied - please adhere to local regulations.						
Surgeons under their own responsibility can hydrate the product with antibiotic solution. The above serves as a reference only, based on published clinical data.						

DAC[®] evidence and references

evidence

HYDROGEL COATING VERSUS CALCIUM SULPHATE BEADS AS A LOCAL ANTIBIOTIC CARRIER FOR DEBRIDEMENT PROCEDURES IN ACUTE PERIPROSTHETIC JOINT INFECTION: A PRELIMINARY STUDY.
Daniele De Meo and Others – MDPI – September 2023

THE INTRAOPERATIVE USE OF DEFENSIVE ANTIBACTERIAL COATING (DAC®) IN THE FORM OF A GEL TO PREVENT PERI-IMPLANT INFECTIONS IN ORTHOPAEDIC SURGERY: A CLINICAL NARRATIVE REVIEW
Daniele Pressato and Others – Materials July 2023

ANTIBIOTIC-LOADED COATINGS TO REDUCE FRACTURE-RELATED INFECTIONS: RETROSPECTIVE CASE SERIES OF PATIENTS WITH INCREASED INFECTIOUS RISK
Daniele De Meo and Others – Antibiotics January 2023

AVOID THE POSSIBILITY OF REMOVE OF OSTEOSYNTHESIS DEVICE SECONDARY TO INFECTION OF THE LUMBAR REGION
Inaki Arroategui, American Journal of Pathology & Research, October 2022

HIGH RATE OF INFECTION ERADICATION FOLLOWING CEMENTLESS ONE-STAGE REVISION HIP ARTHROPLASTY WITH AN ANTIBACTERIAL HYDROGEL COATING
Antonio Pellegrini and Claudio Legnani
The International Journal of Artificial Organs, January 2021
<https://doi.org/10.1177/0391398821995>

CLINICAL APPLICATION OF ANTIBACTERIAL HYDROGEL AND COATING IN ORTHOPAEDIC AND TRAUMATOLOGY SURGERY
Daniele De Meo; Giancarlo Ceccarelli; Giancarlo Iaiani; Federico Lo Torto; Diego Ribuffo; Pietro Persiani, and Ciro Villani
Gels, August 2021
<https://doi.org/10.3390/gels7030126>

ANTIBIOTIC-LOADED HYDROGEL COATING FOR THE PREVENTION OF LOCAL INFECTION AFTER VERTEBRAL SURGERY: A RETROSPECTIVE COHORT ANALYSIS
G. Parbonetti, A. Puglisi, B. Rizzo, E. La Madaia, R. Granata. Surgical Technology International. Volume 39 – August 2021.

RECONSTRUCTION OF INFECTED POST-TRAUMATIC BONE DEFECTS OF THE DISTAL FEMUR WITH THE COMPRESS® IMPLANT. PRELIMINARY RESULTS OF A STAGED NON-BIOLOGICAL STRATEGY.
P. S. Corona, M. Altayó, C. Amat, M. Vicente, R. Velez. Injury. Available online 6 October 2020
<https://doi.org/10.1016/j.injury.2020.10.016>.

ANTIBACTERIAL HYDROGEL COATING IN JOINT MEGA-PROSTHESIS: RESULTS OF A COMPARATIVE SERIES
Carmine Zoccali, Guido Scoccianti, Roberto Biagini, Primo Andrea Daolio, Fabio Luca Giardina, Domenico Andrea Campanacci. The European Journal of Orthopaedic Surgery and Traumatology. Published online February 5th, 2021. <https://doi.org/10.1007/s00590-021-02884-7>

THE POTENTIAL INNOVATIVE USE OF BACTERIOPHAGES WITHIN THE DAC® HYDROGEL TO TREAT PATIENTS WITH KNEE MEGAPROSTHESIS INFECTION REQUIRING “DEBRIDEMENT ANTIBIOTICS AND IMPLANT RETENTION” AND SOFT TISSUE COVERAGE AS SALVAGE THERAPY.
T. Ferry, C. Batailler, C. Petitjean, J. Chateau, C. Fevre, E. Forestier, S. Brosset, G. Leboucher, C. Kolenda, F. Laurent, and S. Lustig. Front. Med. 7:342. July 2020 – doi: 10.3389/fmed.2020.00342

DEFENSIVE ANTIBACTERIAL COATING IN REVISION TOTAL HIP ARTHROPLASTY: NEW CONCEPT AND EARLY EXPERIENCE
M. Franceschini, N.A. Sandiford, V. Cerbone, L. Cappelli Toledo de Araujo, and D. Kendoff. HIP International 2020, Vol. 30(1S) 7–11. DOI: 10.1177/1120700020917125

ANTIBIOTIC-LOADED HYDROGEL COATING TO REDUCE EARLY POSTSURGICAL INFECTIONS IN ASEPTIC HIP REVISION SURGERY: A RETROSPECTIVE, MATCHED CASE-CONTROL STUDY
D. De Meo, V. Calogero, L. Are, A. U. Cavallo, P. Persiani, C. Villani. Microorganisms 2020, 8, 571; doi:10.3390/microorganisms8040571

DEFENSIVE ANTIBACTERIAL COATING (DAC®) FOR PREVENTION OF INFECTION IN ACL RECONSTRUCTION: A FEASIBILITY STUDY
R. Aicale, F. Oliva, N. Maffulli, Muscle, Ligaments & Tendons Journal. Nr 2020;10 (1):151-153. DOI:10.32098/mltj.01.2020.22

ANTIBACTERIAL HYDROGEL COATING IN THE PREVENTION OF PJI AFTER BONE

RECONSTRUCTION WITH MEGAPROSTHESIS: A CONSECUTIVE CASES SERIES; Ep poster accepted at the 2019 EMSOS meeting. Florence 15-17 Maggio 2019.

ANTIBACTERIAL COATING OF IMPLANTS: ARE WE MISSING SOMETHING?
C. L. Romano, H. Tsuchiya, I. Morelli, A. G. Battaglia, L. Drago; Bone Joint Res 2019;8:199–206. doi: 10.1302/2046-3758.85.BJR-2018-0316

TWO-STAGE CEMENTLESS HIP REVISION FOR PERI-PROSTHETIC INFECTION WITH AN ANTIBACTERIAL HYDROGEL COATING: RESULTS OF A COMPARATIVE SERIES.
International Orthopedics – Published online 30 October 2018 – <https://doi.org/10.1007/s00264-018-4206-2>

HYALURONIC BASED ANTIBACTERIAL HYDROGEL COATING FOR IMPLANTABLE BIOMATERIALS IN ORTHOPEDICS AND TRAUMA: FROM BASIC RESEARCH TO CLINICAL APPLICATIONS.
Chapter of the book “Hydrogels” edited by Sajjad Hayder, IntechOpen. Published – August 1st 2018

ONE-STAGE EXCHANGE WITH ANTIBACTERIAL HYDROGEL COATED IMPLANTS PROVIDES SIMILAR RESULTS TO TWO-STAGE REVISION, WITHOUT THE COATING, FOR THE TREATMENT OF PERI-PROSTHETIC INFECTION
Knee Surgery, Sports Traumatology, Arthroscopy. <https://doi.org/10.1007/s00167-018-4896-4>. Published online March 2018

DOES ONE-STAGE EXCHANGE WITH ANTIBACTERIAL COATING OF IMPLANTS PROVIDE SIMILAR RESULTS TO A TWO-STAGE PROCEDURE FOR THE TREATMENT OF PERI-PROSTHETIC JOINT INFECTION?
JBJS VOLUME 99-B, ISSUE SUPP. 1 / JANUARY 2017
Proceedings of the European Orthopaedic Research Society (EORS) 2016, 24th Annual Meeting

ECONOMIC EVALUATION OF ANTIBACTERIAL COATINGS ON HEALTHCARE COSTS IN FIRST YEAR FOLLOWING TOTAL JOINT ARTHROPLASTY
The Journal of Arthroplasty Accepted for publication Feb. 2018

FAST-RESORBABLE ANTIBIOTIC-LOADED HYDROGEL COATING TO REDUCE POST-SURGICAL INFECTION AFTER INTERNAL OSTEOSYNTHESIS: A MULTICENTER RANDOMIZED CONTROLLED TRIAL
J Orthop Traumatol. 2017 Jun;18(2):159-169. doi: 10.1007/s10195-017-0442-2. Epub 2017 Feb 2. ABSTRACT

HYALURONIC ACID AND ITS COMPOSITES AS A LOCAL ANTIMICROBIAL/ANTIADHESIVE BARRIER
J Bone Joint Infect 2017; 2(1):63-72. doi:10.7150/ijbji.17705

HYALURONIC ACID-BASED HYDROGEL COATING DOES NOT AFFECT BONE APPPOSITION AT THE IMPLANT SURFACE IN A RABBIT MODEL
Clin Orthop Relat Res. 2017 Mar 16. doi: 10.1007/s11999-017-5310-0

DOES ONE-STAGE EXCHANGE WITHANTIBACTERIAL COATING OF IMPLANTS PROVIDE SIMILAR RESULTS TO A TWO-STAGE PROCEDURE FOR THE TREATMENT OF PERI-PROSTHETIC JOINT INFECTION?
BJJ ORTHOPAEDIC PROCEEDINGS 2017 vol. 99-B no. SUPP 1 19

ANTI ADHESIVE PROPERTIES OF A FAST-REABSORBABLE, HYDROGEL ANTI-BACTERIAL IMPLANT COATING
BJJ ORTHOPAEDIC PROCEEDINGS 2017 vol. 99-B no. SUPP 2 90

PRINCIPLES OF ORTHOPEDIC INFECTION MANAGEMENT – CHAPTER 9. INFECTION AFTER JOINT ARTHROPLASTY
AOTrauma 2017

DOES AN ANTIBIOTIC-LOADED HYDROGEL COATING REDUCE EARLY POST-SURGICAL INFECTION AFTER JOINT ARTHROPLASTY?
Journal of Bone and Joint Infection – 2016; 1: 34-41. doi: 10.7150/ijbji.15986

ANTIBIOTIC-LOADED HYDROGEL COATING TO PREVENT EARLY POST-SURGICAL INFECTION AFTER JOINT ARTHROPLASTY: RESULTS FROM A MULTICENTRE EUROPEAN TRIAL
BJJ ORTHOPAEDIC PROCEEDINGS 2016 vol. 98-B no. SUPP 23 80

IS ANTIBIOTIC-LOADED HYDROGEL COATING ABLE TO REDUCE EARLY POST-SURGICAL INFECTION AFTER INTERNAL OSTEOSYNTHESIS?
BJJ ORTHOPAEDIC PROCEEDINGS 2016 vol. 98-B no. SUPP 23 43

SYSTEMIC AND LOCAL ADMINISTRATION OF ANTIMICROBIAL AND CELL THERAPIES TO PREVENT METHICILLIN-RESISTANT STAPHYLOCOCCUS EPIDERMIDIS-INDUCED FEMORAL NONUNIONS IN A RAT MODEL
Mediators of Inflammation – Volume 2016 (2016), Article ID 9595706

ANTIBACTERIAL HYDROGEL IMPLANT COATING MAY BE SAFE FOR WRIST, ANKLE OSTEOSYNTHESIS
Orthopaedics Today Europe, May 2015

LOCAL GENTAMICIN DELIVERY FROM RESORBABLE VISCOUS HYDROGELS IS THERAPEUTICALLY EFFECTIVE
Clin Orthop Relat Res. 2015 Jan; 473(1): 337–347. Published online 2014 Sep 17

ANTI-BACTERIAL HYDROGEL COATING OF OSTEOSYNTHESIS IMPLANTS. EARLY CLINICAL RESULTS FROM A MULTI-CENTER PROSPECTIVE TRIAL
European Cells and Materials Vol. 30. Suppl. 2, 2015

THE MANAGEMENT OF PERIPROSTHETIC INFECTIONS IN THE FUTURE: A REVIEW OF NEW FORMS OF TREATMENT
Bone Joint J. 2015 Sep;97-B(9):1162-9

THE MANAGEMENT OF PERIPROSTHETIC INFECTIONS IN THE FUTURE
The Bone & Joint Journal 2015;97-B:1162–9

ANTIBACTERIAL REABSORBABLE COATING OF ORTHOPAEDIC IMPLANTS: AN IN VITRO AND IN VIVO STUDY
BJJ ORTHOPAEDIC PROCEEDINGS 2014 vol. 96-B no. SUPP 11 284

DOES IMPLANT COATING WITH ANTIBACTERIAL-LOADED HYDROGEL REDUCE BACTERIAL COLONIZATION AND BIOFILM FORMATION IN VITRO?
Clin Orthop Relat Res. 2014 Nov; 472(11): 3311–3323. Published online 2014 Mar 13

EFFICACY OF ANTIBACTERIAL-LOADED COATING IN AN IN VIVO MODEL OF ACUTELY HIGHLY CONTAMINATED IMPLANT
Int Orthop. 2014 Jul; 38(7): 1505–1512. Published online 2013 Dec 22

IS VANCOMYCIN-LOADED DAC HYDROGEL COATING OF ORTHOPAEDIC IMPLANTS SAFE FOR HUMAN USE? SHORT-TERM CLINICAL RESULTS IN TWO EUROPEAN CENTERS
Abstract – EBJS 2014

ORTHOPEDIC IMPLANT-RELATED INFECTION PROPHYLAXIS USING A VANCOMYCIN-LOADED, RESORBABLE HYDROGEL COATING. AN IN VIVO ANIMAL STUDY
Abstract – EBJS 2014

NEW HYALURONAN-BASED HYDROGEL AS DISPOSABLE ANTIBACTERIAL COATING FOR PREVENTION OF IMPLANT-RELATED INFECTION IN ORTHOPAEDICS
SIB Meeting 2014

EFFICACY OF ANTIBACTERIAL-LOADED COATING IN ACUTELY HIGHLY CONTAMINATED IMPLANT
SIB meeting 2014

ABSTRACT NO.38370: IN VIVO SAFETY AND EFFICACY OF ANTIBACTERIAL-LOADED DAC HYDROGEL COATING OF ORTHOPAEDIC IMPLANTS
SICOT 2014, Rio De Janeiro, Brazil

RESORBABLE HYDROGEL PROVIDES EFFECTIVE ANTI-BACTERIAL COATING OF IMPLANTS IN VITRO AND IN VIVO
Musculoskeletal Infection Society 2013 Meeting

MEDICATED HYDROGELS OF HYALURONIC ACID DERIVATIVES FOR USE IN ORTHOPEDIC FIELD
Int J Pharm. 2013 Jun 5;449(1-2):84-94

PROPHYLAXIS OF ORTHOPEDIC IMPLANT-RELATED INFECTIONS WITH LOCALLY APPLIED VANCOMYCIN USING A HYDROGEL AS MATRIX
NBTE Lunteren 2013

RESORBABLE HYDROGEL PROVIDES EFFECTIVE ANTI-BACTERIAL COATING OF IMPLANTS IN VITRO AND IN VIVO
Musculo-Skeletal Infection Society annual meeting 2013

NEW HYALURONIC DERIVATE AS ANTIBACTERIAL COATING IN ORTHOPAEDIC SURGERY
Poster presented to the EPNOE 2013 Metting, Nice, France

ABSTRACT NO.35754: DISPOSABLE ANTI-BACTERIAL COATING FOR PREVENTION OF IMPLANT-RELATED INFECTIONS IN ORTHOPAEDICS
Hyderabad OWC 2013

S13.9 ANTIBIOTIC-LOADED REABSORBABLE HYDROGEL COATING FOR INFECTION PROPHYLAXIS OF ORTHOPAEDIC IMPLANTS. PRELIMINARY STUDIES
BJJ ORTHOPAEDIC PROCEEDINGS 2011 vol. 93-B no. SUPP III 337-338

BIOMATERIALS FOR MEDICINE

references

- Kunutsor SK, Whitehouse MR, Blom AW, Beswick AD; Patient-Related Risk Factors for Periprosthetic Joint Infection after Total Joint Arthroplasty: A Systematic Review and Meta- Analysis. PLoS One. 2016;11(3):e0150866. doi: 10.1371/journal.pone.0150866.
- Lenguerrand E, Whitehouse MR, Beswick AD, Kunutsor SK, Burston B, Porter M, Blom AW. Risk factors associated with revision for prosthetic joint infection after hip replacement: a prospective observational cohort study. Lancet Infect Dis. 2018;18(9):1004-1014. doi: 10.1016/ S1473-3099(18)30345-1.
- Assefa M, Amare A. Biofilm-Associated Multi-Drug Resistance in Hospital-Acquired Infections: A Review. Infect Drug Resist. 2022;15:5061-5068. doi: 10.2147/IDR.S379502.
- Sender R, Fuchs S, Milo R. Revised Estimates for the Number of Human and Bacteria Cells in the Body. PLoS Biol. 2016 Aug 19;14(8):e1002533. doi: 10.1371/journal.pbio.1002533.
- Sauer K, Stoodley P, Goeres DM, Hall-Stoodley L, Burmølle M, Stewart PS, Bjarnsholt T. The biofilm life cycle: expanding the conceptual model of biofilm formation. Nat Rev Microbiol. 2022; 20(10):608-620. doi: 10.1038/s41579-022-00767-0
- Pressato D, Battista A, Govoni M, Vivarelli L, Dallari D, Pellegrini A. The Intraoperative Use of Defensive Antibacterial Coating (DAC®) in the Form of a Gel to Prevent Peri-Implant Infections in Orthopaedic Surgery: A Clinical Narrative Review. Materials (Basel). 2023;16(15):5304. doi: 10.3390/ma16155304.
- Parbonetti G, Puglisi A, Rizzo B, La Maida E, Granata R. Antibiotic-loaded Hydrogen Coating for the prevention of local infection after vertebral surgery: A retrospective cohort analysis. Surg Technol Int. 2021;39:441-446. doi: 10.52198/21.STI.39.NS149.
- DAC is a product designed for the preparation and application of a bio-resorbable hydrogel coating mixed with antibiotics. Australian Government TGA - Novagenit Australia Pty Ltd - kit DAC - Implantable-device infection control barrier (427290) - 2023
- Zamboni F, Okoroafor C, Ryan MP, Pembroke JT, Strozyk M, Culebras M, Collins MN. On the bacteriostatic activity of hyaluronic acid composite films, Carbohydr Polym. 2021;260:117803. doi:10.1016/j.carbpol.2021.11780
- Giammona G, Pittarresi G. Palungo SF, Scarponi S , Romano CL. (2018) Hyaluronic-Based Antibacterial Hydrogel Coating for Implantable Biomaterials in Orthopedics and Trauma: From Basic Research to Clinical Applications.Hydrogels. InTech.2018Available at: <http://dx.doi.org/10.5772/>
- Drago L, Boot W, Dimas K, Malizos K, Hänsch GM, Stuyck J, Gawlitta D, Romano CL. Does implant coating with antibacterial-loaded hydrogel reduce bacterial colonization and biofilm formation in vitro? Clin Orthop Relat Res. 2014;472(11):3311-23. doi: 10.1007/s11999-014-3558-1.
- Franceschini M, Sandiford NA, Vincenzo C, Cappelli Toledo de Araujo L, Kendoff D. Defensive antibacterial coating in revision total hip arthroplasty: new concept and early experience. Hip Int. 2020;30(1_suppl):7-11. doi: 10.1177/1120700020917125.
- De Meo D, Calogero V, Are L, Cavallo AU, Persiani P, Villani C. Antibiotic-Loaded Hydrogel Coating to Reduce Early Postsurgical Infections in Aseptic Hip Revision Surgery: A Retrospective, Matched Case-Control Study. Microorganisms. 2020;8(4):571. doi: 10.3390/ microorganisms8040571.
- Trentinagli MT, Van der Streaten C, Morelli I, Logoluso N, Drago L, Romano CL. Economic Evaluation of Antibacterial Coatings on Healthcare Costs in First Year Following Total Joint Arthroplasty. J Arthroplasty. 2018 Jun;33(6):1656-1662. doi: 10.1016/j.arth.2018.01.057.
- Lenguerrand E, Whitehouse MR, Beswick AD, Kunutsor SK, Burston B, Porter M, Blom AW. Risk factors associated with revision for prosthetic joint infection after hip replacement: a prospective observational cohort study. Lancet Infect Dis. 2018;18(9):1004-1014. doi: 10.1016/ S1473-3099(18)30345-1.
- Onggo JR, Onggo JD, de Steiger R, Hau R. Greater risks of complications, infections, and revisions in the obese versus non-obese total hip arthroplasty population of 2,190,824 patients: a meta-analysis and systematic review. Osteoarthritis Cartilage. 2020;28(1):31-44. doi: 10.1016/j.joca.2019.10.00.
- Drago L, Boot W, Dimas K, Malizos K, Hänsch GM, Stuyck J, Gawlitta D, Romano CL. Does implant coating with antibacterial-loaded hydrogel reduce bacterial colonization and biofilm formation in vitro? Clin Orthop Relat Res. 2014;472(11):3311-23. doi: 10.1007/s11999-014-3558-1.
- Pellegrini A, Legnani C. High rate of infection eradication following cementless one-stage revision hip arthroplasty with an antibacterial hydrogel coating. Int J Artif Organs. 2022 Jan;45(1):113-117. doi: 10.1177/0391398821995507.



Distributors



ADLER ORTHO UK
The Stables
Tarvin Road
Frodsham - Cheshire -WA6 6XN
Tel: +44151 329 3372
info@adlerorthouk.com

ADLER ORTHO SpA
Nucleo Produttivo -
Uffici Direzionali
Via dell'Innovazione 9
20032 Cormano (Mi)
Tel +39 02 6154371
Fax +39 02 615437222

ADLER ORTHO SpA
Ufficio di Bologna
Via Guelfa 9
40138 Bologna
Tel. +39-051-533266
Fax +39-051-538772

ADLER ORTHO SpA
Ufficio di Verona
Via Pacinotti 4B
37135 Verona
Tel. +39-045-8230327
Fax +39-045-8230660

ADLER ORTHO SpA
Ufficio di Roma
Via delle Benedettine 86
00135 Roma
Tel. +39 06 30600445
Fax +39 06 30601361

To find out more about how DAC[®] can benefit your practice, please feel free to contact us or visit www.dac-coating.com. All content herein is protected by copyright, trademarks, and other intellectual property rights, as applicable, owned by or licensed to the company or its affiliates unless otherwise indicated. Redistribution, duplication, or disclosure, in whole or in part, is prohibited without the express written consent of the company. This material is intended for healthcare professionals only. The company does not practice medicine. Each physician should exercise independent judgment in the diagnosis and treatment of individual patients. This information does not replace the comprehensive training physicians have received. Check for country-specific regulations and refer to product-specific instructions for use within the packaging for indications, contraindications, warnings, precautions, and potential adverse effects.